

# Validation of a simple RP-HPLC method developed for the quantification of *meta*-cresol in parathyroid hormones formulation

## Abstract

**Aim :** To develop and validate a rapid and sensitive reverse phase high performance liquid chromatography (RP-HPLC) method with UV detection for quantification of *meta*-cresol (*m*-cresol) in pharmaceutical preparation of parathyroid hormone (1–34) (PTH). **Materials and Methods:** Chromatography was performed on a Jupiter RP C-18 (4.6 mm ID × 250 mm L, porosity 300 Å, particle size 5 μm) with a guard column (reversed-phase C<sub>18</sub> column of 4.6 mm ID × 12.5 mm L, porosity 300 Å, particle size 5 μm) using a mobile phase containing 0.1% TFA in 60% methanol with isocratic program at 1.0 mL/min flow rate. Detection was carried out at 217 nm. The method was validated as per ICH guidelines for linearity (correlation coefficient = 0.99), range, accuracy, precision, and robustness ( $n = 9$  during accuracy parameter whereas  $n = 15$  during linearity and range parameter and  $n = 6$  during repeatability). Robustness was confirmed by considering two factors; age effect of the mobile phase and test sample and with different columns during method development. **Results :** The method was linear over the concentration range of 75–120 μg/mL. The precision of the method in terms of relative standard deviation was evaluated from intra- and inter-day replicate injections of system suitability standards of *m*-cresol using different equipment and different columns. Components of within- and between-batch variances were found to be below 2% ( $n = 30$ ) and 3%, respectively, which constituted an acceptable level of variation. Retention time was found to be about 5.2 min and 10.9 min for *m*-cresol and PTH, respectively. **Conclusion:** The developed method thus has the potential of being useful for routine quality control of *m*-cresol.

Key words: *Meta*-cresol, parathyroid hormone, reverse phase HPLC, validation

Shaligram S. Rane,  
Alkesh Ajameri<sup>1</sup>, Rustom  
Mody<sup>1</sup>, P. Padmaja<sup>2</sup>

Departments of Applied  
Chemistry, Faculty of Technology  
and Engineering and Chemistry,  
Science Faculty, The Maharaja  
Sayajirao University of  
Baroda, Vadodara, <sup>1</sup>Intas  
Biopharmaceuticals Ltd.,  
Ahmedabad, Gujarat, India  
<sup>2</sup>Department of Chemistry,  
Science Faculty, The Maharaja  
Sayajirao University of Baroda,  
Vadodara 39002, Gujarat, India.  
E-mail: p\_padmaja2001@  
yahoo.com

### Address for correspondence:

Dr. Padmaja Sudhakar,  
Department of Chemistry,  
Science Faculty, The Maharaja  
Sayajirao University of Baroda,  
Vadodara 39002, Gujarat, India.  
E-mail: p\_padmaja2001@yahoo.  
com

### Access this article online

Website: [www.phmethods.org](http://www.phmethods.org)

DOI: 10.4103/2229-4708.90366

Quick response code



## INTRODUCTION

*Meta*-cresol (*m*-cresol) is used as bactericide in the biotechnological processing of pharmaceuticals; preservative in pharmaceutical formulations [injection solutions of insulin, somatotropin, and parathyroid hormone (PTH)]; pesticide for the treatment of the stems of fruit trees and plants. Exposure of humans is possible through the use of *m*-cresol as a preservative in pharmaceutical injection solutions.

*Meta*-cresol, *para*-cresol, and *m/p*-cresol mixtures are absorbed across the respiratory and gastrointestinal tracts and through the skin, and are distributed throughout the body. The primary metabolic pathway for all cresol isomers is conjugation with glucuronic acid and inorganic sulfates. All isomers are mainly eliminated by renal excretion in the form of above-mentioned conjugates. The oral LD<sub>50</sub> of undiluted *m*-cresol in rats was 242 mg/kg bw. Clinical sign includes hypoactivity, salivation, tremors, and convulsions. Neither mortality nor clinical signs of toxicity were seen following exposure to saturated vapor concentration of either *m*-cresol or *p*-cresol. Inhalation of aerosols may, however, cause death, and mean lethal concentrations in rats were reported to be 29 mg/m<sup>3</sup> for *p*-cresol and 58 mg/m<sup>3</sup> for *m*-cresol.<sup>[1]</sup> Reaction to *m*-cresol in commercial preparations of insulin to humans was reported by Dennis *et al*<sup>[2]</sup>.

The analysis of cresol-like chemicals in use for a long period of time has evolved from a number of nonspecific colorimetric methods to more selective separation techniques using gas chromatography (GC) or high performance liquid chromatography (HPLC).<sup>[3-5]</sup>

The objective of our work was hence to develop a rapid and simple RP-HPLC method with UV detection, useful for routine quality control of *m*-cresol in PTH formulations. Human PTH, a peptide of 84 amino acid residues<sup>[6]</sup> secreted from parathyroid gland, is the principle homeostatic regulator of the level of blood calcium through its actions on kidney and bone.<sup>[7]</sup> Teriparatide (recombinant DNA origin) injection [recombinant human PTH (1–34)] is a bone-forming agent used for the treatment of osteoporosis. The structure of meta-cresol and PTH are given in Figure 1 (a) and (b), respectively.

The method developed was validated for parameters such as linearity, accuracy, and precision. So far, to the best of our knowledge, no RP-HPLC method has been reported for determination of *m*-cresol in PTH formulations.

## MATERIALS AND METHODS

### Material, reagent, and chemicals

HPLC-grade acetonitrile and methanol were purchased from Merck; tri-fluoro-acetic acid was purchased from Sigma Aldrich. Ultra pure water was obtained using Milli-Q® UF-Plus (Millipore) system; *m*-cresol was obtained from J.T. Baxter/Hedinger and was used for preparation of different dilutions; PTH API (Active Pharmaceutical Ingredient) having concentration of 400 µg/mL was used for different dilutions of a PTH working standard. PTH

formulation containing 250 µg/mL PTH as active pharmaceutical ingredient was used as a test sample. All chemicals, i.e. mannitol, sodium acetate, and glacial acetic, were of the highest purity available.

### Preparation of standard, mobile phase, and dilution buffer

#### Formulation buffer

Buffer containing 3 mg/mL *m*-cresol, 45.4 mg/mL mannitol, 0.1 mg/mL sodium acetate and 0.41 mg/mL glacial acetic acid in Milli Q water was prepared. It is similar as excipients used for PTH formulation.

#### PTH working standard

PTH (400 µg/mL) was used for preparation of different diluted samples.

#### *m*-cresol standard

An aliquot of 3 mg/mL was used for preparation of different dilutions.

#### Mobile phase consisted of 0.1% v/v TFA in 60% methanol

All dilutions were made using calibrated digital micro-pipettes.

### Chromatographic conditions

LC system equipped with an injection valve (quaternary), 217 UV detector, and Chemstation software was used. A reversed-phase Jupiter C<sub>18</sub> column (4.6 mm ID × 250 mm L, porosity 300 Å, particle size 5 µm) with a guard column (reversed-phase C<sub>18</sub> column of 4.6 mm ID × 12.5 mm L, porosity 300 Å, particle size 5 µm) was used for separation. To get the optimum results, the mobile phase with a flow rate of 1.0 mL/min was used and the column temperature was maintained at 30°C. The isocratic programme for the mobile phase was optimized for 12 min.

### Validation of chromatographic methods

During validation, an analytical method is confirmed for reliability, accuracy, and precision of the intended purpose of that method.<sup>[8]</sup> The International Conference on Harmonisation (ICH) of Technical Requirements for Registration of Pharmaceuticals for Human Use<sup>[9]</sup> recognizes accuracy, precision, repeatability, intermediate precision, specificity, limit of quantitation and detection, linearity, and the range of a method as important validation parameters. The guidance for the methodology and statistical tests and associated limits for the validation procedure is provided in the ICH Q2B.<sup>[10]</sup>

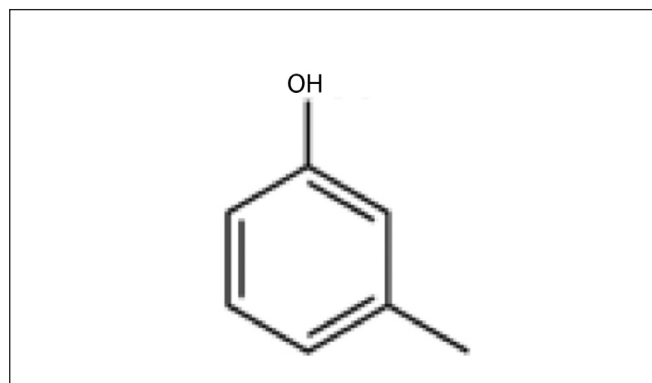


Figure 1: Structure of *m*-cresol

**Table 1: Linearity and range results**

Sample name	Conc. (mg/mL)	Principal peak area			Average total area	% RSD	Back calculated conc. (mg/mL)	% recovery
		Prep. 1	Prep. 2	Prep. 3				
Standard-1	0.075	2460	2366	2410	2412	1.9	0.04	100.33
Standard-2	0.091	2910	2899	2899	2903	0.2	0.05	98.53
Standard-3	0.100	3230	3354	3315	3300	1.9	0.07	101.33
Standard-4	0.110	3534	3657	3596	2596	1.7	0.10	100.02
Standard-5	0.120	3841	3972	3963	2925	1.9	0.50	99.75

**Table 2: Results of different test parameters**


Statistical parameter	Details/results
Linearity and range (conc. in µg/mL)	75–120
Regression equation	$y = 34039x - 149.44$
Correlation coefficient	0.9969
Total analysis time in min	12 min
Retention time in min	
For <i>m</i> -cresol	About 5.2
For PTH	About 10.9
Flow rate	1 mL/min
Column	C18 (4.6 × 250 mm) 5 µm, 300 Å, Guard column
Column condition	Flow rate = 1.0 mL/min
Mobile phase	0.1% TFA in 60% methanol
Method	Isocratic
Validation parameter	
Specificity	No interference
Accuracy	Recovery more than 95%
Precision	RSD less than 2%
Robustness	Method is robust for all tested parameters

The optimized chromatographic methods were validated according to the procedures described in ICH guidelines Q2 (R1).<sup>[11]</sup> The USP<sup>[12]</sup> states that the validation of an analytical method is a process that establishes the performance characteristics of a developed analytical method and ensures that it meets its intended purpose and analytical application. Method validation includes an assessment of the adequacy of the analytical procedure by means of statistical testing, including linear regression analysis, and relative standard deviation (RSD) determination in order to demonstrate the validity of the method.<sup>[13]</sup>

## RESULTS AND DISCUSSION

### Method development

To obtain the best chromatographic conditions, the mobile phase composition, column temperature, and flow rate were adequately selected. The flow rate was varied from 0.8 mL/min to 1.2 mL/min. The column temperature was varied between 22°C and 30°C, and the analysis at 30°C was preferred to improve

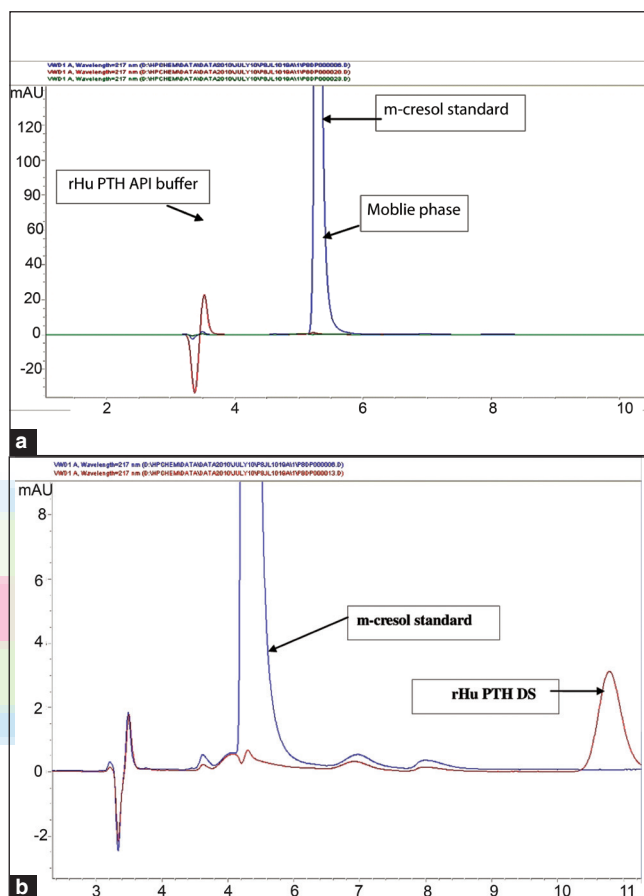

**Figure 2:** (a) Sample chromatogram of the principal peak of PTH API buffer, *m*-cresol standard and mobile phase. (b) sample chromatogram (overlap) of the principal peak of PTH API and *m*-cresol standard

the peak symmetry and resolution. The % of the mobile phase was varied from 58% to 62%. Isocratic chromatographic conditions were optimized for determination of *m*-cresol in a PTH pharmaceutical product. The chromatographic separation was achieved by applying chromatographic conditions described in “Chromatographic condition” section.

The applied chromatographic conditions permitted a good separation of *m*-cresol and PTH at different concentrations of *m*-cresol. No interference of other excipients was observed as shown in Figure 2 a and b.

The capacity factor ( $k'$ ) of the first peak (*m*-cresol) and second peak (PTH) was 3.24 and 5.24, respectively;

**Table 3: Accuracy results**

Sample name	Conc. (mg/mL)	Area of principal peak area (mAu s)				% RSD	Back calculated conc. (mg/mL)	% recovery
		Prep. 1	Prep. 2	Prep. 3	Average			
Standard-1	0.091	3032	3028	3025	3028	0.1	0.093	102
andard-2	0.100	3183	3182	3179	3181	0.1	0.098	98
andard-3	0.120	4030	4031	4027	4029	0.1	0.121	101


**Table 4: Repeatability results**

Area of principal peak of six independent prep. (mAu s)	Average area of principal peak (mAu s)	RSD of principal peak of six injections (%)	Rt of principal peak (min)
3154	3262	1.8	5.3
3241			5.3
3302			5.3
3257			5.3
3308			5.3
3310			5.3

while the resolution factor was 6.88. The asymmetry of the peak for *m*-cresol and PTH were found to be 1.29 and 5.29, respectively; while the tailing factor parameter for *m*-cresol and PTH was found to be 1.29 and 1.14, respectively. For replicate injections of *m*-cresol standard; the % RSD of the main peak area was found to be below 0.7%, and there was no significant variation in the retention time (<0.1 min).

The PTH and *m*-cresol peaks were thus found to be well resolved, and the tailing factor was within the limits. Available PTH formulation in the market contains 100 µg/mL of *m*-cresol and the range was selected as 75–120 µg/mL. Detection limit (LOD) and quantification limit (LOQ) are not required for the study.

Different brands of reverse phase C<sub>18</sub> columns were used (Jupiter column and Grace Vydac column) and compared in terms of percentage variation of principal peak area of *m*-cresol standard. Experiments were conducted using a system suitability standard. Percentage variation between principal peak of *m*-cresol was not more than 5% in all samples when compare to the area of principal peak of *m*-cresol from specificity samples. Retention time of the principal peak of *m*-cresol was found to be around 5.3 min and 3.9 min whereas the principal peak of PTH was found to be around 10.9 and 4.7 min on Jupiter and Grace Vydac columns, respectively. Relative retention rime with Grace Vydac and Jupiter columns of PTH with respect to *m*-cresol was found to be 1.2 and 2.1, respectively. The principal peak in both samples was separated by base-to-base while overlapping their chromatograms. On the basis of relative retention time, a Jupiter column was selected for method validation.

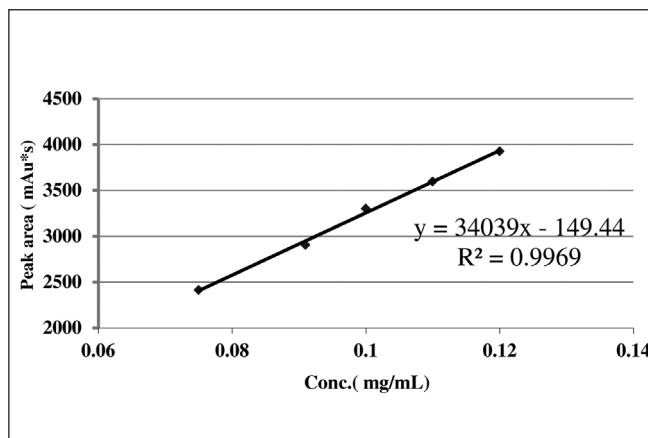**Figure 3: Linearity**

## Method validation

### Specificity

Specificity of the method for *m*-cresol in the presence of PTH, and excipients such as mannitol, sodium acetate, and glacial acetic acid was studied in terms of resolution of peaks observed and peak area of *m*-cresol. To verify this, PTH, formulation buffer, mobile phase, and *m*-cresol standard were injected onto HPLC separately. Three different concentrations of *m*-cresol (75, 100, and 120 µg/mL) were prepared in the mobile phase as well as in the formulation buffer and were injected in triplicate onto HPLC. From Figures 2(a and b) it is evident that there is no interference and the method developed is specific for *m*-cresol.

### Linearity and range

The ICH guidelines for an assay method recommend that 80–120% of stated value can be considered as a range of a method. We decided to study linearity in the concentration range from 75 to 120 µg/mL as the concentration of *m*-cresol in the drug product is 100 µg/mL. *m*-Cresol standard (3 mg/mL) was used for preparation of different concentrations ranging from 75 to 120 µg/mL, by considering 100 µg/mL as 100%. Five different concentrations were considered with three replicates of each concentration (n = 15). A linearity curve was plotted for peak area responses versus concentration of *m*-cresol as shown in Figure 3, and the results are shown in Table 1. The correlation coefficient, slopes, Y-intercepts, and the regression equation of the calibration curve were determined and shown in

**Table 5: Experimental design for intermediate precision**

	Set 1	Set 2	Set 3	Set 4	Set 5
Equipment used	HPLC equipment-1	HPLC equipment-2	HPLC equipment -2	HPLC equipment-1	HPLC equipment-1
Day	Day-1	Day-2	Day-2	Day-3	Day-3
Column	Column-1	Column-2	Column-1	Column-1	Column-1

Table 2. The % RSD was found to be <2.0% while the % recovery was found to be in the range of 97–103%.

#### Accuracy

Accuracy was studied using two different sets of three different solutions, containing 90, 100, and 120 µg/mL of *m*-cresol. Each solution was spiked in the mobile phase and injected onto HPLC (n = 9). The % recovery was found to be between 98% and 102% and the % RSD was found to be <1.0% as shown in Table 3.

#### Precision

Precision was evaluated based on intra-day (repeatability) and inter-day (intermediate precision) variation and on different columns. The repeatability was assessed with six independent samples of 100 µg/mL of *m*-cresol. Single injection from each preparation was injected and the results are shown in Table 4. The % RSD of the main peak area was found to be <2.0 %.

Using the experimental design and matrix as shown in Table 5, intermediate precision was evaluated on different days with different equipment and with different columns. Three replicate injections of system suitability standards (100 µg/mL of *m*-cresol) prepared independently were considered for the study. Intra-day precision was determined for 100 µg/mL of *m*-cresol by performing five different conditions as mentioned in Table 5 (n = 15) and RSD were calculated.

The % RSD for the main peak area of *m*-cresol standard within each set and between different sets was found to be <2.0%. The % recovery of *m*-cresol standard was found to be between 95.0% and 105.0% within each set, and the maximum variation between sets was found to be 3.0%. The results are shown in Tables 1, 3, and 4.

#### Robustness

Robustness is a measure of its capacity to remain unaffected by small, but deliberate variations in method parameters and provides an indication of its robustness during normal usage. Robustness was tested using two variables: age effect of the mobile phase and test samples.

Freshly prepared samples (100 µg/mL of *m*-cresol) and those stored for 7 days were analyzed using both

the freshly prepared mobile phase and mobile phase stored for 7 days. No variation in retention time was observed with percent variation from the initial day to 7 days. However, percentage variation of the principal peak area of *m*-cresol was found to be higher than 5% when compared to the principal peak area of *m*-cresol in specificity samples. Due to variation in the peak area, it was decided to use freshly prepared sample as well as mobile phase before analysis.

## CONCLUSION

The RP-HPLC method was demonstrated to be validated for quantifying *m*-cresol in the presence of other excipients. Validation studies showed that the developed HPLC method was selective, linear, precise, and accurate. The chromatographic method described here was found to be reliable for quantifying *m*-cresol in PTH formulation. Since the method is simple and rapid, it may be successfully applied to quality control analysis of *m*-cresol in PTH formulations.

## REFERENCES

1. m/p Cresol category, UNEP Publication. Available from: <http://www.inchem.org/documents/sids/sids/m-p-cresols.pdf> [Last accessed on 2011 Aug 10].
2. Dennis K. Baraniuk J. Delayed-type hypersensitivity reaction to the *m*-cresol component of insulin. *Ann Allergy Asthma Immunol* 2007;99:194-5.
3. McNeil D. Kirk-Othmer Encyclopedia of Chemical Technology". Vol. 6, 2<sup>nd</sup> ed. New York: John Wiley and Sons; 1965. p. 434-44.
4. Husain S, Kunzelmann P, Schildknecht H. Separation of isomeric alkyl phenols by high-performance liquid chromatographic and gas-liquid chromatographic techniques. *J Chromatogr* 1977;13:753-60.
5. "NIOSH Manual of Analytical Methods". Vol. 3, 2<sup>nd</sup> ed. USDHEW, PHS, CDC, NIOSH, DHEW (NIOSH) Publication No. 77 - 157C, April 1977.
6. Hendy GN, Kronenberg HM, Potts JT Jr, Rich A. Nucleotide sequence of cloned cDNAs encoding human preproparathyroid hormone'. *Proc Natl Acad Sci USA* 1981;78:7365-9.
7. Potts JT Jr, Kronenberg HM, Rosenblatt M. Rosenblatt, Parathyroid hormone: chemistry, biosynthesis, and mode of action. *Adv Protein Chem* 1982;35:323-96.
8. US Food and Drug Administration, Center for Drug Evaluation and Research, Reviewer Guidance, Validation of Chromatographic Methods, Rockville, MD, 1994. <http://www.fda.gov/CDER/GUIDANCE/cmc3.pdf>. [Last accessed on 2005 Apr 05].
9. Code Q2A - Text on Validation of Analytical Procedures International Conference on Harmonization (ICH) of Technical Requirements for the registration of Pharmaceuticals for Human Use, Guideline for Industry; 1995. Available from: <http://www.fda.gov/cder/Guidance/ichq2a.pdf>. [Last accessed on 2005 Apr 05].

**Rane, et al.: RP-HPLC method for quantification of *m*-cresol**

10. Code Q2B - Validation of Analytical Procedures: Methodology, International Conference on Harmonization (ICH) of Technical Requirements for the registration of Pharmaceuticals for Human Use, Guidance for Industry; 1996. Available from: <http://www.fda.gov/cder/Guidance/1320fnl.pdf>. [Last accessed on 2005 Apr 05].
11. Code Q2-R1 - Validation of analytical procedures: Text and methodology, ICH harmonized tripartite guideline; 2005.
12. United States Pharmacopeia, 29<sup>th</sup> ed. United States Pharmacopeial Convention, Inc, Twin brook Parkway, Rockville, MD, USA; 2006. p. 1623, 2675-93, 3050-3, 3392-4.
13. US Food and Drug Administration, Center for Drug Evaluation and Research, Guidance for Industry: Analytical Procedures and Methods Validation. Chemistry, Manufacturing, and Controls

Documentation DRAFT GUIDANCE, Rockville, MD, 2000. Available form: <http://www.fda.gov/CDER/GUIDANCE/2396dft.pdf>. [Last accessed on 2005 Apr 05].

**How to cite this article:** Rane SS, Ajameri A, Mody R, Padmaja P. Validation of a simple RP-HPLC method developed for the quantification of *meta*-cresol in parathyroid hormones formulation. Pharm Methods 2011;2:203-8.

**Source of Support:** NMPB, Government of India (Project No. GO/MH-04/2009), **Conflict of Interest:** None declared.

